BRAUN

ThermoScan™ 7


Ear thermometer

IRT6520B

GB	English	2
DE	Deutsch	9
ES	Español	16
FR	Français	23
IT	Italiano	30

NL	Nederlands	37
PL	Polski	44
SA	اللغة العربية	57
SE	Svenska	58

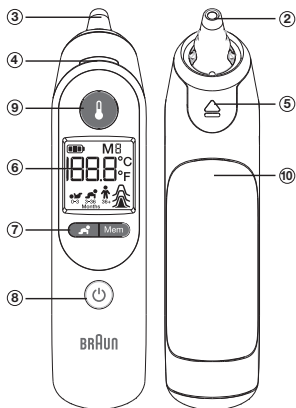Certain trademarks used under license from The Procter & Gamble Company or its affiliates.

ThermoScan™, ExacTemp™, and Age Precision™ are trademarks of Kaz Europe Sàrl and/or its affiliates.

Kaz Europe Sàrl
A Helen of Troy Company
Q-Center - Chaux 4
CH-1030 Bussigny - Switzerland

Made in Vietnam
Printed in China

© 2026, All rights reserved. www.BraunHealthcare.com/uk_en/

Patents: www.BraunHealthcare.com/us_en/patents

Made in Vietnam

MD REF IRT6520BWEV1

P/N: A013886R0
22OCT25

Product description

- ① Disposable Hygiene cap (a box of 20 is supplied + 1 already installed) LF40/LF20
- ② Probe tip
- ③ Probe
- ④ ExacTemp™ Stability Indicator
- ⑤ Hygiene cap ejector
- ⑥ Display
- ⑦ Age Precision™ Technology and Memory buttons
- ⑧ Power button
- ⑨ Start button
- ⑩ Battery door (2 X 1.5V type AA (LR06) batteries)
- ⑪ Protective case

Product description

The Braun ThermoScan™ Ear thermometer has been carefully developed for accurate, safe, and fast temperature measurements in the ear. The shape of the thermometer probe prevents it from being inserted too far into the ear canal which can hurt the eardrum. However, as with any thermometer, proper technique is critical to obtaining accurate temperatures. Therefore, read the instructions carefully and thoroughly.

Intended use and contraindications

The Braun ThermoScan™ Ear thermometer is intended to be used by home users for intermittent measurement of the body temperature of patients having ages ranging from normal weight (full term) new-born to geriatric adult. The thermometer is reusable. The Braun Hygiene cap is a single use device and is intended to maintain accuracy of measurement, hygiene and prevent cross-contamination.

This product is not intended to diagnose any disease but is a useful screening tool for temperature. Use of this thermometer is not intended as a substitute for consultation with your doctor.

This thermometer is not intended for clinical use in a professional environment and is for home usage only. The Age Precision™ Technology feature is not intended for pre-term babies or small-for-gestational age babies. This thermometer is not intended to interpret hypothermic temperatures. Do not allow children under 12 to take their temperatures unattended.

Use of the IRT6520 is not intended as a substitute for consultation with your physician.

Safety information and important precautions

- Please consult your physician if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc. in spite of any colour of the background light or absence of fever.
- Please consult your physician if the thermometer shows elevated temperature (yellow or red background lights).

The operating ambient temperature range for this thermometer is 10 – 40 °C (50 – 104 °F). Do not expose the thermometer to temperature extremes (below -25 °C / -13 °F or over 55 °C / 131 °F) or excessive humidity (>95% RH). This thermometer must only be used with genuine Braun ThermoScan™ disposable Hygiene caps (LF 40).

To avoid inaccurate measurements always use this thermometer with a new, clean Hygiene cap attached.

If the thermometer is accidentally used without a Hygiene cap attached, clean the lens (see «Care and cleaning» section). Keep Hygiene caps out of reach of children.

Do not modify this equipment without authorization of the manufacturer.

Parents/guardians should call the pediatrician upon noticing any unusual sign(s) or symptom(s).

For example, a child who exhibits irritability, vomiting, diarrhea, dehydration, seizure, changes in appetite or activity, even in the absence of fever, or who exhibits a low temperature, may still need to receive medical attention.

Children who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

Temperature elevation as indicated by Age Precision™ may signal a serious illness, especially in adults who are old, frail, have a weakened immune system, or neonates and infants. Please seek professional advice immediately when there is a temperature elevation and if you are taking temperature on:

- neonates and infants under 3 months (consult your doctor immediately if the temperature exceeds 37.4 °C or 99.4 °F)
- patients over 60 years of age
- patients having diabetes mellitus or a weakened immune system (e.g. HIV positive, cancer chemotherapy, chronic steroid treatment, splenectomy)
- patients who are bedridden (e.g. nursing home patient, stroke, chronic illness, recovering from surgery)
- a transplant patient (e.g. liver, heart, lung, kidney)

Fever may be blunted or even absent in elderly patients.

This thermometer contains small parts that can be swallowed or produce a choking hazard to children. Always keep the thermometer out of children's reach.

Age Precision™ Technology

Clinical research shows the definition of fever changes as newborns grow into children and children grow into adults¹.

Braun's Age Precision™ Technology takes the guesswork out of interpreting temperature for the entire family.

Simply use the Age Precision™ Technology button to select the appropriate age setting, take a temperature and the colour-coded display will glow green, yellow or red to help you understand the temperature reading.

¹ Herzog L, Phillips S. Addressing Concerns About Fever. *Clinical Pediatrics*. 2011; 50(#5): 383-390.

How to use your Braun ThermoScan™

1. Remove the thermometer from the protective case.

2. Push the Power button .


During an internal self-check, the display shows all segments. Then the last temperature taken will be displayed for 5 seconds.

3. The Hygiene cap indicator will blink to signal one is needed.

To achieve accurate measurements, make sure a new, clean Hygiene cap is in place before each measurement.

Attach a new Hygiene cap by pushing the thermometer probe straight into the Hygiene cap inside the box and then pulling out.


Note: The Braun ThermoScan™ will not work unless a Hygiene cap is attached.

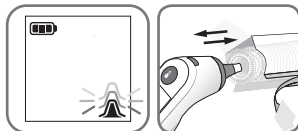4. Select age with Age Precision™ Technology button  as shown below. Press button to toggle through each age group.

NOTE: You must select an age in order to take a measurement.


0-3 Months

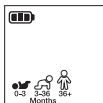
3-36 Months

36 Months – adult

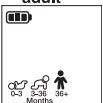
The thermometer is ready to take a temperature when the display looks like the images below.

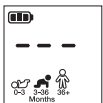5. Fully insert the IRT probe into the ear canal and point it towards the ear canal, make sure the IRT is pointed toward the ear drum and not the canal walls.

Once positioned correctly, press and release the start button .

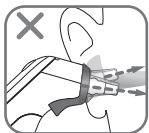6. ExacTemp™ Stability Indicator light will pulse while temperature measurement is in progress. The light will remain solid for 3 seconds to indicate that a successful temperature reading has been achieved.




NOTE: If the probe has been properly inserted into the ear canal during the measurement, a long beep will sound to signal the completed measurement.

If the probe has not been constantly placed in a stable position in the ear canal, a sequence of short beeps will sound, the ExacTemp™ Stability Indicator light will go out and the display will show an error message (POS = position error).

See «Errors and troubleshooting» section for more information.

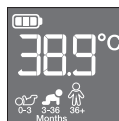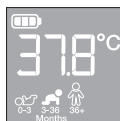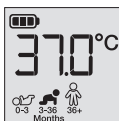7. The confirmation beep indicates that an accurate temperature measurement has been taken. The result is shown on the display. Temperature reading displayed is an equivalent of oral temperature.


Age range	Green Normal temperature	Yellow Fever	Red High fever
 0-3 Months	≥35.8 – ≤37.4 °C (≥96.4 – ≤99.4 °F)	X	>37.4 °C (>99.4 °F)
 3-36 Months	≥35.4 – ≤37.6 °C (≥95.7 – ≤99.6 °F)	>37.6 – ≤38.5 °C (>99.6 – ≤101.3 °F)	>38.5 °C (>101.3 °F)
 36 Months – adult	≥35.4 – ≤37.7 °C (≥95.7 – ≤99.9 °F)	>37.7 – ≤39.4 °C (>99.9 – ≤103.0 °F)	>39.4 °C (>103.0 °F)

8. For the next measurement, press Eject button  to remove and discard used Hygiene cap, and put on a new, clean Hygiene cap.

NOTE: Defaults to last age setting used if you do not change.

The Braun ThermoScan™ Ear thermometer turns off automatically after 60 seconds of inactivity. The thermometer can also be turned off by pressing the Power button .

The display will briefly flash OFF and it will go blank.

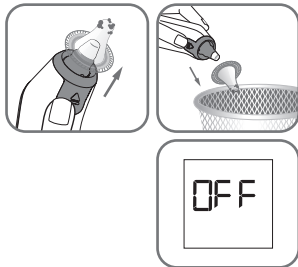Temperature-taking hints


Replace disposable Hygiene cap with every use to maintain accuracy, hygiene and to prevent cross-contamination. Earwax and dirt on the tip of the Hygiene cap, even in a microscopic layer, can create an obstacle between the sensor and the infrared heat emitted by the ear and impact the accuracy of the measurement.

External factors like poor probe replacement, used Hygiene cap or dirty Hygiene cap may influence ear temperature measurements.

In clinical studies, precision and accuracy of ear thermometer measurement has been shown to be influenced by variability in user technique. **Correct probe placement is essential for accurate measurements.**

Temperature in right and left ear differs, always compare measurements from the same ear. The natural physiological difference may reach 0.2 °C. In the absence of confounding factors, in high-risk population (like infants or immunodepressed people) use the highest of the two to interpret the patient's condition.

In the following cases, wait 20 minutes prior to taking a temperature: extreme changes in room temperature, hearing aid, lying on a pillow. Use the untreated ear if ear drops or other medications have been placed in the ear canal.

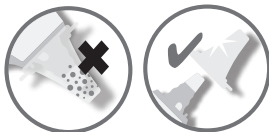
Memory mode

1. This model stores the last 9 temperature measurements. To display the stored measurement, the thermometer must be turned on.


2. Press the **Mem** button.

The display shows the temperature, and when releasing the **Mem** button, the stored temperature for that memory number is displayed, together with "M".

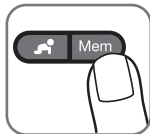3. Each time the **Mem** button is pushed, the stored reading and an M is displayed to indicate each temperature measurement (eg. M2).


Memory mode is automatically exited by not pressing the memory button for 5 seconds.

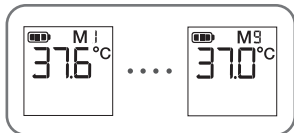How to illuminate the display

The thermometer includes a convenient night light to illuminate the display in a dim environment.

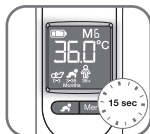
The light will turn on when you press any button. It will stay on until the thermometer is inactive for up to 15 seconds even after temperature is taken.

How to change temperature between °C and °F

Your Braun ThermoScan™ is shipped with the Celsius (°C) temperature scale activated. If you wish to switch to Fahrenheit (°F) and/or back from Fahrenheit to Celsius, proceed as follows:

1. Make sure the thermometer is turned off.
2. Press and hold down the Power button . After about 3 seconds the display will show this sequence: °C / SET / °F / SET.
3. Release the Power  button when the desired temperature scale is shown. There will be a short beep to confirm the new setting, then the thermometer is turned off automatically.

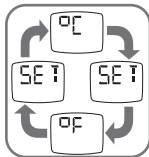Care and cleaning

 The probe tip is the most delicate part of the thermometer. It must be clean and intact to ensure accurate readings. If the thermometer is ever accidentally used without a Hygiene cap, clean the probe tip as follows:

Very gently wipe the surface with a cotton swab or soft cloth moistened with alcohol. After the alcohol has completely dried out, you can put a new Hygiene cap on and take a temperature measurement.

If the probe tip is damaged, please contact Consumer Relations.

Use a soft, dry cloth to clean the thermometer display and exterior. Do not use abrasive cleaners. Never submerge this thermometer in water or any other liquid. Store thermometer and Hygiene caps in a dry location free from dust and contamination and away from direct sunlight.

Additional Hygiene caps (LF 40) are available at most stores carrying Braun ThermoScan™.

Replacing the batteries

The thermometer is supplied with two 1.5 V type AA (LR 06) batteries. For best performance, we recommend alkaline batteries.

1. Insert new batteries when the battery symbol appears on the display.
2. Open the battery compartment. Remove the batteries and replace with new batteries, making sure the poles are in the right direction.
3. Snap battery door into place.

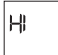
 Only discard empty batteries. They should not be disposed of in the household waste, but at appropriate collection sites or at your retailer.

Errors and troubleshooting

Error message	Situation	Solution
	No Hygiene cap is attached.	Attach new, clean Hygiene cap.
	The thermometer probe was not positioned securely in the ear. An accurate measurement was not possible.	Take care that the positioning of the probe is correct and remains stable. Replace Hygiene cap and reposition. Press Start button to begin a new measurement.

	Ambient temperature is not within the allowed operating range (10 – 40 °C or 50 – 104 °F).	Allow the thermometer to remain for 30 minutes in a room where the temperature is between 10 and 40 °C or 50 and 104 °F.
	Temperature taken is not within typical human temperature range (34 – 42.2 °C or 93.2 – 108 °F).	Make sure the probe tip and lens are clean and a new, clean Hygiene cap is attached. Make sure the thermometer is properly inserted. Then, take a new temperature.
	HI = too high LO = too low	
	System error – self-check display flashes continuously and will not be followed by the ready beep and the ready symbol. If error persists, If error still persists,	Wait 1 minute until the thermometer turns off automatically, then turn on again. ... reset the thermometer by removing the batteries and putting them back in. ... please contact the Service Centre.
	Battery is low, but thermometer will still operate correctly.	Insert new batteries.
	Battery is too low to take correct temperature measurement.	Insert new batteries.

Product specifications

Displayed temperature range:	34 – 42.2 °C	(93.2 – 108 °F)
Operating ambient temperature range:	10 – 40 °C	(50 – 104 °F)
Storage temperature range:	-25 – 55 °C	(-13 – 131 °F)
Operating relative humidity:	10-95% RH (non condensing)	
Storage relative humidity:	15-95% RH	
Display resolution:	0.1 °C or °F	

Accuracy for displayed temperature range

35 – 42 °C (95 – 107.6 °F):	Maximum Laboratory Error	±0.2 °C	(±0.4 °F)
Outside this range:		±0.3 °C	(±0.5 °F)
Clinical repeatability:		±0.14 °C	(±0.26 °F)
Battery life:	2 years / 350 measurements		
Battery Type:	2 X 1.5V type AA (LR 06) batteries		
Service life:	5 years		
Accessory:	Hygiene cap, LF40/LF20, single use		
Reference Body Site:	Oral Equivalent		

This thermometer is specified to operate at 1 atmospheric pressure or at altitudes with an atmospheric pressure up to 1 atmospheric pressure (700 - 1060 hPa).

Refer to Instruction manual (blue symbol)

Equipment with type BF applied parts

Consult electronic instruction for use

Operating temperature

Storage temperature

95%
10%
Operating humidity

95%
15%
Storage humidity

Unique device identifier

Batch code

Serial number

Catalogue number

Medical device

Hygiene cap: Do not re-use

Manufacturer

Keep it dry

Authorized representative in the European Community

Subject to change without notice.

This medical device carries the CE mark and is manufactured in conformity with RoHS Directive 2011/65/EU and other applicable directive and/or regulations as specified in the EU Declaration of Conformity. Full responsibility for the conformance of the product to the standard is assumed by Kaz Europe Sarl, Place Chauderon 18, CH-1003 Lausanne, Switzerland.

ASTM laboratory accuracy requirements in the display range of 37 to 39 °C (98 to 102 °F) for IR thermometers is ± 0.2 °C (± 0.4 °F), whereas for mercury-in-glass and electronic thermometers, the requirement per ASTM Standards E667-86 and E1112-86 is ± 0.1 °C (± 0.2 °F).

EMC Tables are available on our website at www.hot-europe.com/support.

The device needs special precautions regarding EMC. Portable and mobile RF communications equipment can affect medical electrical equipment. **WARNING:** Use of this equipment adjacent to other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally. **WARNING:** Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30cm to any part of the device. Otherwise, degradation of the performance of this equipment could result. **WARNING:** use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at your local retailer or at appropriate collection points provided in your country.

Calibration

This product has been calibrated at the time of manufacture. If used according to the use instructions, periodic calibration is not required. If at any time, you question the measuring accuracy, please contact our service representative (see last page for contact information) or visit www.HelenofTroy.com/emea-en/support.

Inquiry

In Europe, any serious incident (e.g., death, life-threatening injury, surgical intervention, etc.) that has occurred in relation to this thermometer should be reported to Kaz Europe (see Helen of Troy consumer service numbers) and to the competent authority of your member state.

REF, LOT, and Serial numbers of your product must be provided as they are essential to record and follow up your inquiry or claim. The Serial number SN is located in the battery compartment. The production date is given by the LOT number, which is located in the lower housing of the device, and can be deciphered as explained below:

The first 3 numeric digits in the LOT number represent the day of the year of manufacture. The next 2 numeric digits represent the last two numbers of the calendar year of manufacture and the letter(s) at the end designate the manufacturer of the product. (e.g.: LOT No.: 12313tav this product was made on the day 123, year 2013 at the manufacturer code tav.)

Guarantee

For support requests, please visit www.HelenofTroy.com/emea-en/support or find service contact information on the last page of this manual.

UK only: This guarantee in no way affects your rights under statutory law.

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